

M55B/C



Physical Assessment Simulator "Physiko" Plus

*only M55C has the CE certificate

Learning-to-Listen features open wider possibilities for training

"Learning-to-Listen" training: teaching medical students to anticipate, describe and identify physical examination findings.

The goal of "Learning-to-Listen" is to provide learners with deliberate practice in the recognition and interpretation of clinical findings obtained during the auscultation of the heart sound.

Product Supervision

Toyoaki Yamauchi, MD, ND, PhD
Professor
Department of Nursing Graduate School of
Medicine Nagoya University



KEY FEATURES

Learning-to-Listen

Nine modules with contrasting scenarios to facilitate learning for all medical and nursing students.

- In each module, preset patients are prepared for individual physical assessment.
- After the physical assessment, trainees can take programmed heart sound test.

Select a module




Select a patient



Physical assessment



Heart sound test

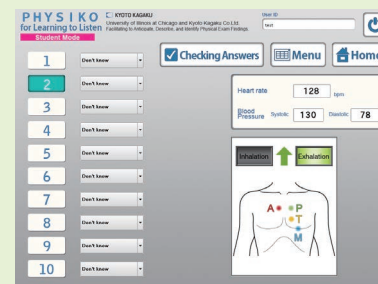


Physical Assessment Mode 1


Patient Cases

Physiko's twelve preset patient scenarios cover typical complaints such as thoracic pain and breathing difficulty.

Top menu

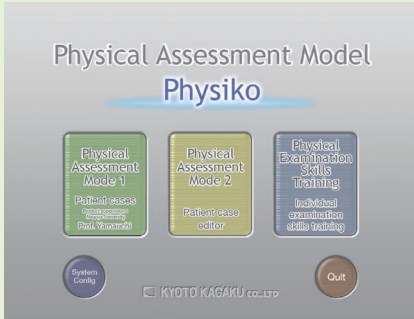

Choose a patient

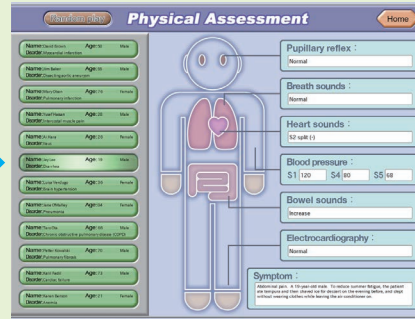


Start a session



Cases Comprehensive examination of 12 typical patients

Case	Age, Sex	Condition
Thoracic Pain	58, Male	Myocardial Infarction
	55, Male	Dissecting Aortic Aneurysm
	76, Female	Lung Infarction
	28, Male	Intercostal Muscle Ache
Abdominal Pain	26, Female	Ileus
	19, Male	Diarrhea


Case	Age, Sex	Condition
Unconscious	36, Female	Brain Hypertension
	84, Female	Pneumonia
	66, Male	Chronic Obstructive Lung Disease (COPD)
Stuffy Feeling	70, Male	Pulmonic Fibrosis
	73, Male	Heart Failure
	21, Female	Anemia

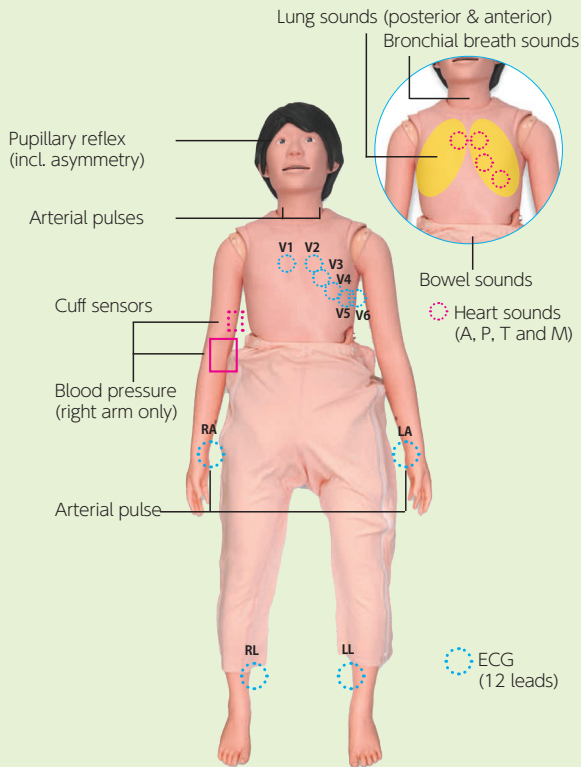

Physical Assessment Mode 2

Patient Cases Editor


Users can create new simulated patients by selecting from the built-in options.

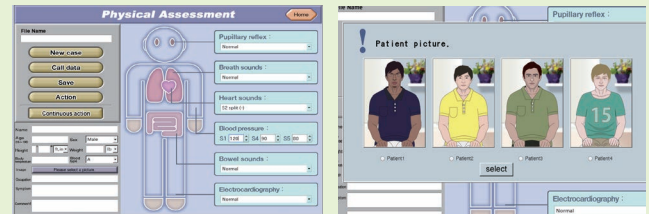

Physical Assessment Skills Training

Individual Examination Skills Training

Physiko allows individual assessment skills training

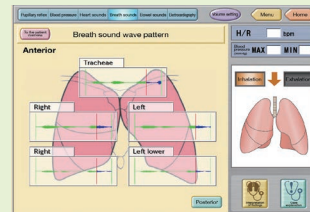

- Pupillary reflex
- Blood pressure
- Auscultation of lung sounds
- Auscultation of heart sounds
- Auscultation of bowel sounds
- ECG simulation

DESCRIPTIONS

SPECIFICATIONS

Size: 155cm/61 inch Weight: 24kg/53lbs
Power: AC100V-240V 50/60Hz
Power Consumption: 70W

SET INCLUDES

- 1 manikin
- 1 terminal box
- 1 ECG electrodes unit*
- 1 laptop PC
- 1 sphygmomanometer unit
- 1 set of clothes
- 1 microphone

*ECG electrodes unit is not included in M55C

MATERIALS

Soft resin, hard resin
Latex free