



Height measuring instruments

for measuring body height

5003

Height rod, foldable

SOEHNLE
PROFESSIONAL

- S** Handy height rod for mobile use
- S** For measuring yourself and others

- > Ultrasound measuring technology, temperature-compensated
- > Tilt sensor for correction if measuring instrument is not horizontal, resulting in an error message
- > Low space requirement due to fold function
- > Plastic housing
- > LCD display, 9 mm digit height
- > Acoustic acknowledgment at end of measurement and reading displayed
- > Units in cm or inch
- > Battery operation with 2 x 1.5 V AA batteries (included in scope of delivery)



| | |
|-----------------|------------------|
| Order number | 5003.01.001 |
| Measuring range | 50 - 240 cm |
| Division | 1 cm |
| Dimensions | 345 x 92 x 28 mm |
| Scale weight | 250 g |



Folded perfect for mobile use.



5002.01

Measuring tape

S *Small and handy*

- > Measuring tape for wall-mounting
- > No wear parts and maintenance-free
- > Units in mm

| | |
|-----------------|---------------|
| Order number | 5002.01.001 |
| Measuring range | 0 - 220 cm |
| Division | 1 mm |
| Dimensions | 150x30x140 mm |
| Scale weight | 130 g |



5005


Baby height rod

SOEHNLE
PROFESSIONAL

S *Digital and analogue readings*

- > Easy to use, quick and precise
- > Hold function: The reading remains on the display
- > LCD display, 16 mm digit height
- > Incl. adapter for attachment to baby scale 7725 (included in scope of delivery)
- > Units in mm or inch
- > Energy-saving auto-off function
- > Battery operation with 2 x 1.5 V AAA batteries (included in scope of delivery)

| | |
|-----------------|---------------|
| Order number | 5005.10.001 |
| Measuring range | 35 - 80,1 cm |
| Division | 1 mm |
| Dimensions | 630x290x70 mm |
| Scale weight | 600 g |

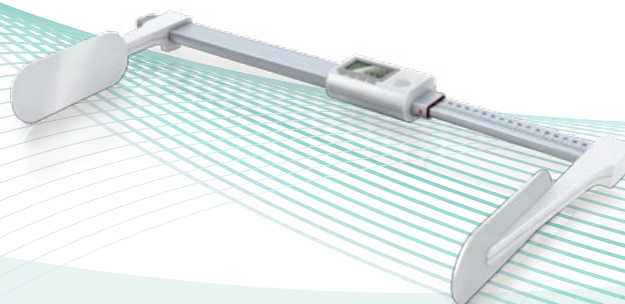

5002.02

Height rod

S *Simple wall-mounting system*

- > Height rod for wall-mounting
- > Headstop hinged
- > No wear parts and maintenance-free
- > Mounting kit included in scope of delivery

| | |
|-----------------|----------------|
| Order number | 5002.02.001 |
| Measuring range | 105 - 210 cm |
| Division | 0.5 cm |
| Dimensions | 50x1,175x22 mm |
| Scale weight | 800 g |

