

Solutions for Hospital-use Blood Pressure Monitors

Mercury-Free
Professional Blood Pressure Monitors

Durability & Easy-to-Clean

The body and cuff of the UM Series are chemical resistant and easy to clean with the antiseptic solutions (such as ethanol, chlorhexidine gluconate and benzalkonium chloride).



Fits Various Patients

The newly designed professional cuff of the UM Series ranges from SS to LL (5 sizes), 12cm to 50cm for the arm circumference.

Cuff Size	Arm Size	REF (for UM-102)	REF (for UM-201& UM-211*)
SS	12cm to 17cm	CUF-KW-SS	CUF-KS-SS
Small Adult	16cm to 24cm	CUF-KW-SA	CUF-KS-SA
Adult	22cm to 32cm	CUF-KW-A	CUF-KS-A
Large Adult	31cm to 45cm	CUF-KW-LA	CUF-KS-LA
LL	41cm to 50cm	CUF-KW-LL	CUF-KS-LL

*Adult, Large Adult and Small Adult (3 sizes) are available for use with UM-201.



Easy to Store & Easy to Carry

In order to safely and securely carry the device, a cuff holder and grip for carrying are available for the UM Series. Highly convenient to carry from room to room, patient to patient.

○ Cuff Holder ○ Grip for Carrying



Specifications

	UM-102	UM-201	UM-211
Product picture			
Measurement method	Auscultatory measurement with stethoscope	Oscillometric measurement	Oscillometric measurement and Auscultatory measurement with stethoscope
Measurement range	Pressure: 0 - 300mmHg Pulse: 40 - 180 beats / minute	Pressure: 0 - 299mmHg Pulse: 40 - 180 beats / minute	Pressure: 0 - 299mmHg Pulse: 40-200 beats / minute
Measurement accuracy	Pressure: ±3mmHg, Pulse: ±5%	Pressure: ±3mmHg, Pulse: ±5%	Pressure: ±3mmHg, Pulse: ±5%
Power supply	2 x 1.5V Batteries (LR6 or AA)	4 x 1.5V Batteries (LR6 or AA)	AC 100V - 240V (Rechargeable Battery)
Battery life	Approx. 2,000 measurements	Approx. 450 measurements	Approx. 300 measurements
Memory storage	N / A	N / A	Last 99 measurements
Operating conditions	+10°C to 40°C / 15%RH to 85%RH/800 hPa to 1060 hPa	+10°C to 40°C / 15%RH to 85%RH/800 hPa to 1060 hPa	+5°C to 40°C / 10%RH to 85%RH/800 hPa to 1060 hPa
Transport / Storage conditions	-20°C to +60°C / 10%RH to 95%RH	-10°C to +60°C / 15%RH to 85%RH	-20°C to +60°C / 10%RH to 95%RH
Dimensions	Approx. 98<W> x 67<H> x 324<D> mm (A-type) Approx. 98<W> x 326<H> x 202<D> mm (B-type)	Approx. 150<W> x 156<H> x 120<D> mm	Approx. 120<W> x 200<H> x 140<D> mm
Weight	Approx. 520g (Unit without Battery) (A-type) Approx. 540g (Unit without Battery) (B-type)	Approx. 480g (Unit without Battery)	Approx. 550g (Unit without Battery)
Standard accessories	Adult Cuff (22-32cm), Rubber Bulb Unit, 2 x AA Battery	Adult Cuff (22-32cm), 4 x AA Batteries	Adult Cuff (22-32cm), Battery Pack, AC Adaptor, AC Cable
Options	CUF-KW-LL: LL Cuff (41-50cm) CUF-KW-LA: Large Adult Cuff (31-45cm) CUF-KW-A: Adult Cuff (22-32cm) CUF-KW-SA: Small Adult Cuff (16-24cm) CUF-KW-SS: SS Cuff (12-17cm) UM-102-10: Rubber Bulb Unit UM-102-11: Connector and Plugs for Cuff Replacement UM-ST001: Mobile Stand (for B-type only)	CUF-KS-LA: Large Adult Cuff (31-45cm) CUF-KS-A: Adult Cuff (22-32cm) CUF-KS-SA: Small Adult Cuff (16-24cm) TB-240: AC Adaptor (A Plug) TB-233C: AC Adaptor (C Plug) TB-233BF: AC Adaptor (BF Plug)	CUF-KS-LL: LL Cuff (41-50cm) CUF-KS-LA: Large Adult Cuff (31-45cm) CUF-KS-A: Adult Cuff (22-32cm) CUF-KS-SA: Small Adult Cuff (16-24cm) CUF-KS-SS: SS Cuff (12-17cm) UM-211-20: Battery Pack TB-268: AC Adaptor K01886: AC Cable (A Plug) K01887: AC Cable (C Plug) K01888: AC Cable (BF Plug)

HPRA Wholesalers authorisation No. W00426/00001
Tel: +353 (0) 1 835 2411 Email: sales@medguard.ie
Unit 28B Ashbourne Business Centre, Ashbourne, Co Meath, A84 WA49

AND
A&D Medical

MERCURY FREE

MEDGUARD
Professional Healthcare Supplies



UM-102

UM-201

UM-211

CE 0123

AND
A&D Medical

Technology You Can Trust

<http://www.aandd.jp>



'Minamata Convention on Mercury'

In October 2013, The 'Minamata Convention on Mercury' was adopted in Minamata city (in Japan where serious health incident occurred as a result of mercury pollution in the mid-20th Century). Under this treaty, the manufacture, import or export of mercury and mercury contained equipment will be prohibited globally in 2020. The 'Minamata Convention on Mercury' requires Mercury Blood Pressure Monitors have to be phased out.

A&D Medical is Responding to 'Mercury-Free Era'

Under the global movement towards Mercury-Free products, replacing Mercury Blood Pressure Monitors with alternative and reliable Mercury-Free devices will be a significant issue for all daily users. As a leading manufacturer of digital blood pressure monitors, A&D is responding to the 'Minamata Convention on Mercury', and has developed a series of new professional blood pressure monitors, the UM Series.

UM Series Lineup



UM-102 A-type
Auscultatory Hybrid
Traditional Style Mercury-Free Sphygmomanometer



UM-102 B-type
Auscultatory Hybrid
Mercury-Free Sphygmomanometer with New Dual Use Design



UM-201
Oscillometric
Simple Operation Blood Pressure Monitor



UM-211
Oscillometric and Auscultatory
All-in-one Blood Pressure Monitor


UM-102

Clinically Validated CE 0123
Mercury-Free Sphygmomanometer

In 2005, A&D developed the world's first hybrid* traditional style sphygmomanometer. UM-102 is the 2nd generation model in this category. Two separate versions of the UM-102 are available. We provide a standard portable model as well as a dual-use model.

*Mercury-Free + Auscultatory Measurement + Digital Display

- ◆ Auscultatory blood pressure measurement
- ◆ Digital readout for pulse rate
- ◆ Cuff holder and grip for carrying
- ◆ Durable and chemical resistant body
- ◆ Full ABS housing light weight body
- ◆ Professional cuff and 5 cuff sizes are available for use
- ◆ Adjustable display column for easier analysis of result
- ◆ Easy battery replacement
- ◆ Anti-glare display minimizes reflected light

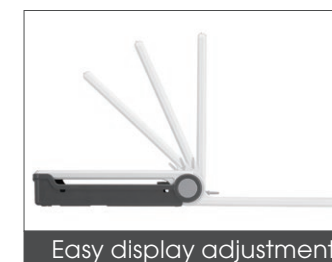

Easy display adjustment



Easy battery replacement

A-type Traditional Style



A-type



Flat base for stacking multiple units



B-type Stand-alone use



In combination with optional mobile stand

B-type New Dual Use Design

UM-201

IHB Clinically Validated CE 0123
Simple Operation Blood Pressure Monitor

Currently many home-use blood pressure monitors are used in hospitals because of their reasonable cost and simple operation. UM-201 is our optimal solution to meet the required durability at hospitals and replace blood pressure monitors designed for home-use. This easy-to-use device comes with a clear display, durable body and professional cuff.

- ◆ Durable and chemical resistant body
- ◆ Professional cuff and 3 cuff sizes are available for use
- ◆ One-touch operation
- ◆ Large LCD display
- ◆ Cuff holder and grip for carrying
- ◆ Irregular HeartBeat (IHB) indicator



UM-211

IHB Clinically Validated CE 0123
All-in-one Blood Pressure Monitor

UM-211 is designed for multipurpose use in hospitals. This model is available for dual measurement mode (auscultatory mode with use of the stethoscope and oscillometric mode).

- ◆ Dual measurement mode (Auscultatory and Oscillometric)
- ◆ Cuff holder and grip for carrying
- ◆ Durable and chemical resistant body
- ◆ Professional cuff and 5 cuff sizes are available for use
- ◆ Large LCD display with backlight
- ◆ Rechargeable battery (300 measurements at full charged)
- ◆ Adjustable inflation pressure (AUTO/220/250/280)
- ◆ Irregular HeartBeat (IHB) indicator
- ◆ Room temperature display

UM-211's auscultatory mode is designed for daily and frequent use. Selectable deflation speed (2.5mmHg / sec. or 5.0mmHg / sec.) and quick start for auscultatory mode are features 'Only from A&D'.



Stand-by display (Room temperature)



Measurement display (with backlight)

Multiple Usages



Outpatient



Stored and charged at nurse station



Easy to carry



Measurement at ward

