


## YXH-4A/B/C Scoop Stretcher

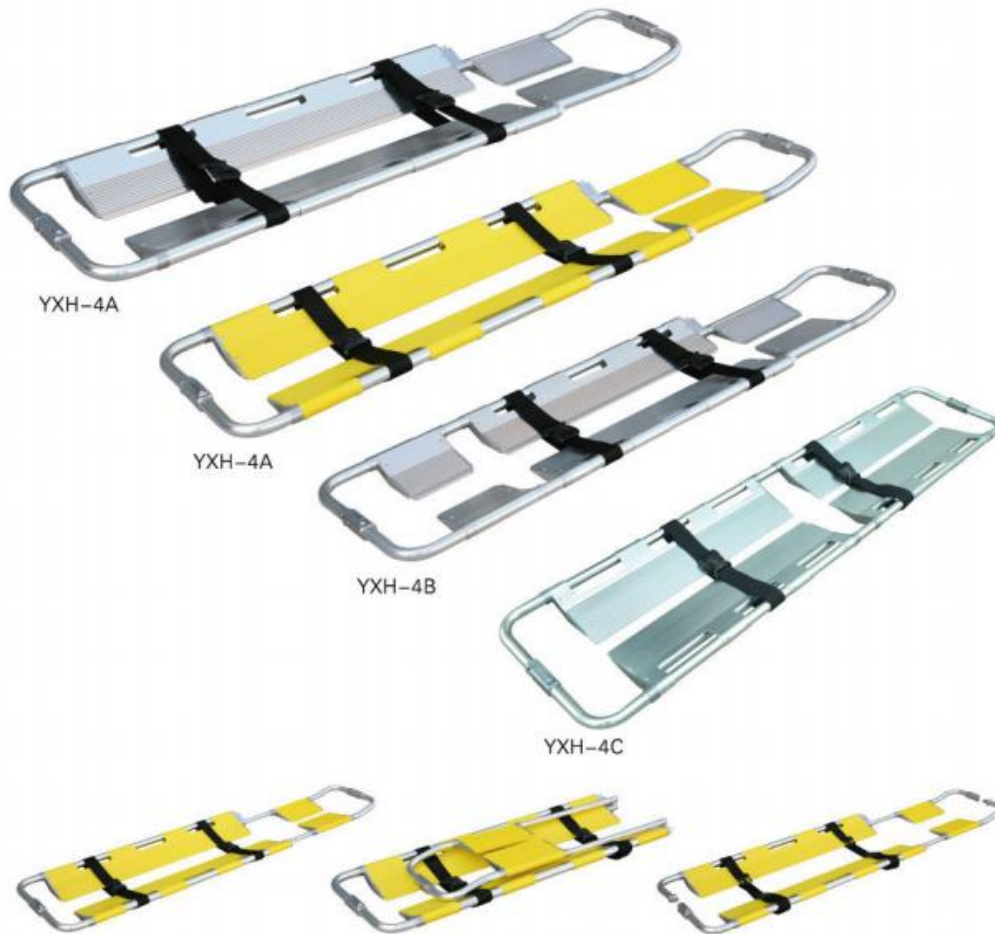

### Detailed Product Description:

It is a separately-type stretcher to transfer fractured and serious patients. At both ends of the stretcher, a hinged and joined arrangement is set in the middle, which can set the stretcher into two parts and it is easy to lock and unlock. The stretcher can strengthen the patients at original place; it can reduce the further hurt of the patients, the patients can be quickly & easily put in the stretcher and also the stretcher can be taken out from the back of the patient without moving the patient. The length of the stretcher can be adjusted according to patients' body. The foot end is designed as a narrow-frame structure. It is made of high-strength aluminum alloy materials. It characterized by its being light-weighted, small-sized, easy-carried, use-safely and easy for sterilization and clean. It is mainly used for hospitals, sports, ambulance carrying patients and wounded person. It also has yellow and blue color.

Model	Product Size (L*W*H)	Folded Size (L*W*H)	Packing Size (1pcs)	Load Bearing	N.W	G.W
YXH-4A/B	210 x 44 x 6cm	120x 44 x 9cm	170 x 45 x 7cm	≤159kg	8.5kg	10.2kg
YXH-4C	185 x 44 x 6cm		188 x 46 x 8cm	≤159kg	9.2kg	11kg