



Overview:

The G3 model features electric operation of height, backrest angle and Trendelenburg tilt facility for improved patient positioning for gynaecology and colposcopy procedures. This couch has a seat section with a cut-out area and the option of a removable extension which increases the overall length of the couch. The Trendelenburg facility provides pelvic tilt to -20° and offers a recovery position in the event of fainting or light-headedness. Additionally, the electric backrest can be positioned up to 85° for sitting procedures or lowered to -15° to aid patient recovery when required. Standard accessories include the matching head cushion, height adjustable armrests, leg supports, and removable debris tray.


Benefits

- Lifetime frame & 5 year parts guarantee - see website for further details
- Electric operation of height, backrest and tilt facility for ease of patient positioning
- Leg supports with both height adjustment and multi-angular placement
- Seat section with cut-out area and option of extension for patient positioning and recovery
- Height adjustable arm supports for patient comfort and security
- Safe working load of 250kgs (550lb / 39 stone)
- Retractable wheels provide easy movement and situating of the couch
- Electrically adjustable backrest angle from $+85^{\circ}$ down to -15° to aid patient recovery
- Low patient transfer height of 52cm (20")
- Wide choice of options, width and accessories to meet user requirements

Overall Dimensions

Length: 144cm (57") approx
with seat extension 180cm (71") approx
Width: 63cm (25") approx
Height: 52cm - 101cm approx
(20" - 40") approx

Upholstery Dimensions

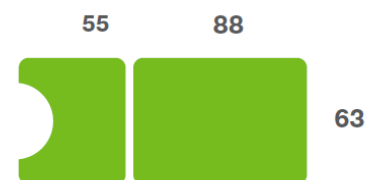


Options & Accessories:




Classic Operators Chair




Foot Switch Operation

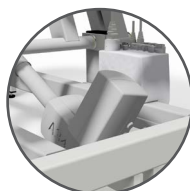


Battery Back-up

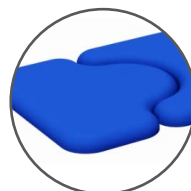

Seat Extension

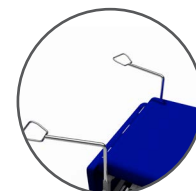


Lithotomy Stirrups

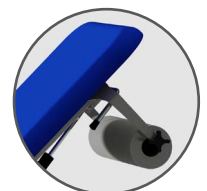

Paper Roll Holder

Operational Parameters

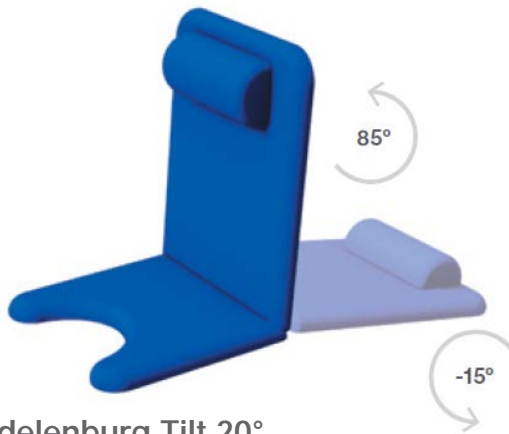

Trendelenburg Tilt 20°



Model 93G3 - Electric Version

Standard Features

- Designed and manufactured in the UK
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Electric height operation with Trendelenburg tilt
- Comfort foam
- Easy operation
- Low running costs
- Hardwearing anti-microbial vinyl
- Safe working load 250kg
- Retractable castors
- Hoist clearance 130mm
- Leg supports
- Arm loops
- Head cushion
- Stainless tray
- Motorised backrest

Technical Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601 P (EN)
- Vinyl Flame Retardant to CRIB 5 (BS 7176)
- Foam to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]
- CE Marked
- Externally audited to ISO9001:2015 by LRQA

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service may be available if required
- Overseas delivery 7-10 days plus shipping time

Upholstery Vinyl

- Flame retardant to BS 5852 Ignition Source 5 - CRIB 5 (BS 7176) Medium Hazard
- Independently tested to verify effectiveness at reducing the spread of MRSA, E.coli & other similar bacteria & fungi
- Rated at 500,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high-traffic, heavy-duty areas
- Phthalate free
- Please note: If the upholstery colour is not specified, 'Almond' will be supplied



Jasmine



Battleship



Cobalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond

Medguard Professional Healthcare

Tel: +353 (0) 1 835 2411

Email: sales@medguard.ie

Fax: 01 969 5009

Website: www.medguard.ie

Unit 6/7, Block 4, Ashbourne Business Park, Ashbourne, Co. Meath, A84 VK52