

LUPA LED H.F

LED magnifier with round lens, Lupa LED H.F is intended when performing any meticulous work in need of great accuracy. It is ideal for hospitals, measured studies, laboratories, aesthetic clinics, etc.

Functionality and Performance

Generates a light without dark spots.
Internal spring balanced arms that allow different positions.
Switch on the head.
Round head with convex lens, adjustable and simple mechanism

Mounting Options

The standard version includes AH table clamp.

As for accessories, we offer a wide selection of supports such as wall supports, wall extension arms, trolleys stand and rail brackets.

Cleaning and Safety

Head exempt of heat.
Protective cover on the head with hinge mechanism to prevent from dust and scratches that may damage the lens.

Economy

High performance, low power 8 W.
Long life, 40,000 hours.
Durable materials.
Very low maintenance costs.



LUPA LED H.F

Technical Details	(REF. 09340)
Light source	LED H.F 8W
Useful life	40.000 hours
Color temperature	6.500 °K
Lenses	5 diopters
Magnification Lens Diameter	127 mm
Articulated Arm Length	900 mm
Lighting unit Net Weight	2,40 kg
Color	White



TROLLEY STAND
8,8 kg REF. 09651






WALL B SUPPORT
REF. B0101800




AH TABLE CLAMP
REF. B010190

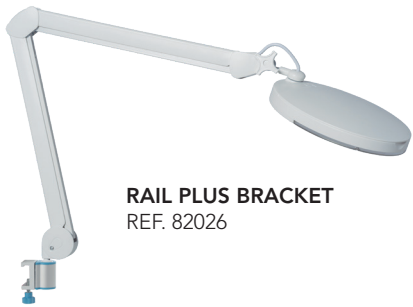

RAIL PLUS BRACKET
REF. 82026




EXTENSION ARM
REF. B0101S0

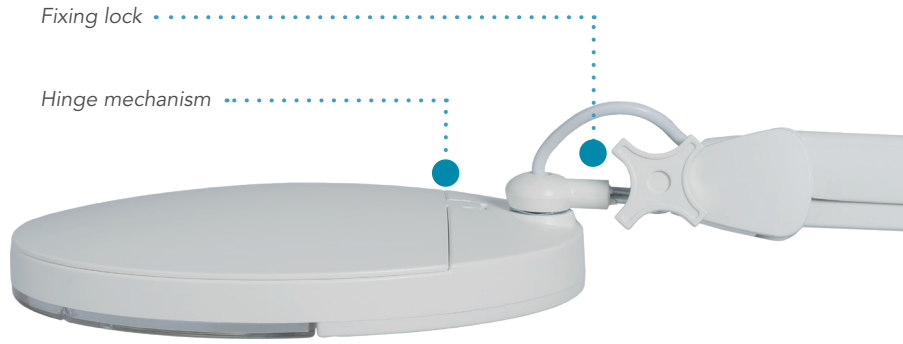


Fixing lock

Hinge mechanism

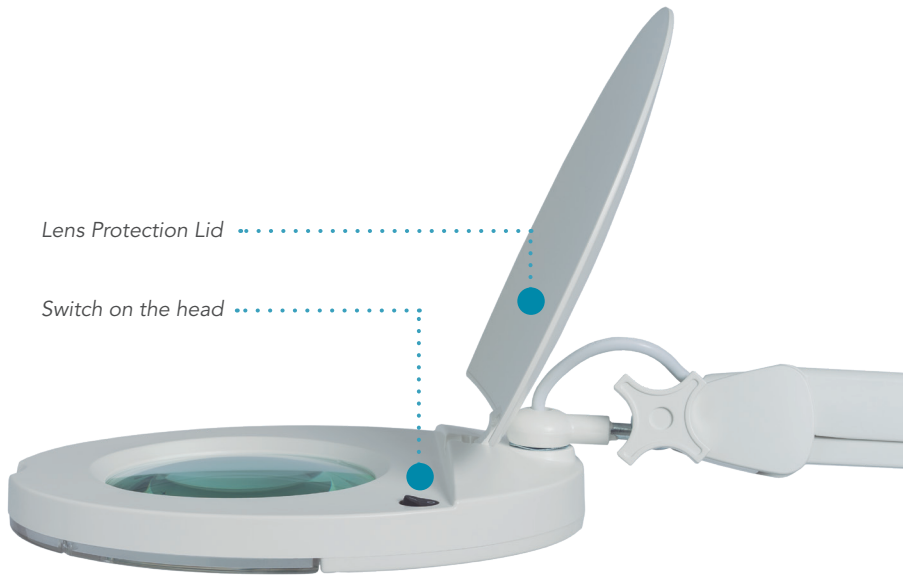

Lens Protection Lid

Switch on the head

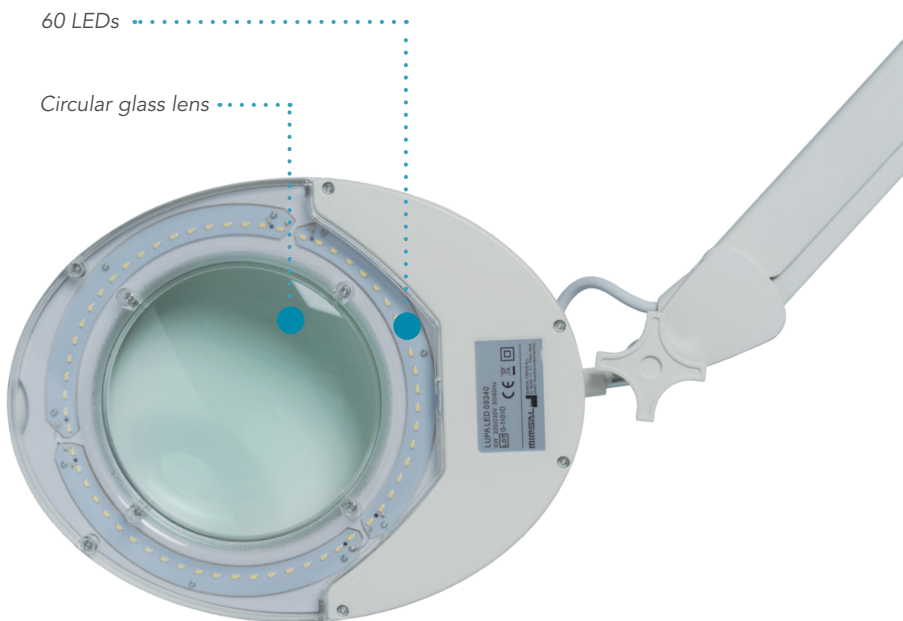

60 LEDs

Circular glass lens

Medguard Professional Healthcare

Tel: 01 835 2378

Email: orders@medguard.ie

Fax: 01 969 5009

Web: www.medguard.ie

Unit 6/7, Block 4, Ashbourne Business Park, Ashbourne, Co. Meath, Ireland