


MW56

Laparoscopic Suture Training System -Laparo Cube



Quantitative Feedback combined with Realistic Hands-on Experience

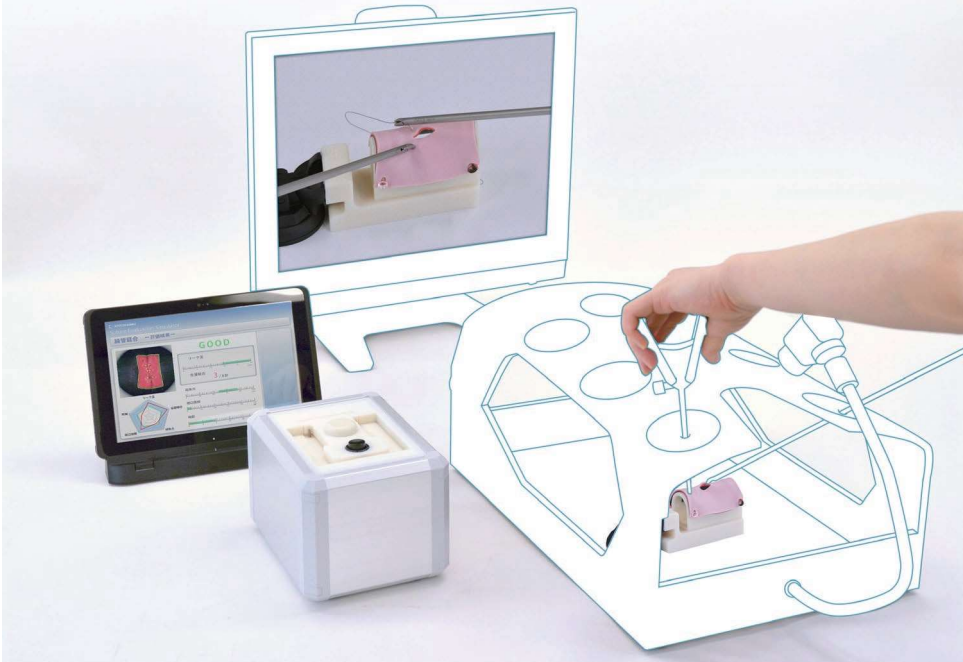

FEATURES

- 1 | Innovative assessment criteria and quantitative measurement provide objective feedback to trainers.
- 2 | Trainees' personal skills assessment data can be saved for reviewing and analysis.
- 3 | Actual clinical instruments can be used.
- 4 | Realistic and lifelike suturing sheet.
- 5 | Able to create your own passing standard or adjust criteria threshold to trainees' level.

SKILLS

- | Laparoscopic suturing (simple interrupted suture)

KEY FEATURES

Judging success or failure

- 1 | Is the wound sutured without a leak?
- 2 | Is suture completed through muscular layer to mucosal layer?
- 3 | Is ligature appropriate for suture?
- 4 | How much open area is there in mucosal layer?
- 5 | How long does it take for suturing procedure?



Life-like Intestinal Tract

Intestinal tract has four layers with respective characteristics. When holding with forceps and suturing, these four layers react like real tissue.

Easy to set up!



DESCRIPTIONS

SPECIFICATIONS

Size: W15 x D20 x H13 cm /
W5.9 x D7.9 x H5.1 inch

Weight: 1.5kg / 3.3lbs
(Assessment unit)

Power: AC100V-240V 50/60Hz
Power Consumption:
26VA : Assessment unit
40VA : PC

MATERIALS

Soft resin
Latex free

RECOMMENDED DEVICES

-

SET INCLUDES

1 laparoscopic suturing assessment unit	1 PC
1 power cable for assessment unit	1 power cable for PC
1 set of skin sheet holder	30 pieces of suturing sheet (intestinal tract)
1 cable	1 container
1 USB Cable	1 instruction manual

REPLACEMENT PARTS

11395-010 30 pieces of suturing sheet
※ 1 set consist of inner layer& outer layer

Box trainer, display monitor, and forceps and camera are not included in set