

seca 675

Electronic platform scales III



- Generously sized platform and very high capacity.
- Suitable for weighing persons in a wheelchair or on a chair.
- Flexible cable remote display.
- Transport castors.

seca 675:

Weighing without restrictions

Special situations demand special scales. Not every patient can step onto a scale quickly and easily. Most of all, patients with limited mobility or in wheelchairs and weak or seriously overweight persons need the help of a convenient scale like the seca 675. The platform scale serves a wide range of patients and eases the routine work of doctors and nurses.

Large platform with high load-bearing capacity.

Simply weigh the wheelchair too: The oversized non-slip platform of the seca 675 offers plenty of space for all conventional wheelchair models. For patients who have trouble standing, a chair can be placed on the scale – comfort and convenience that cannot be weighed. The pre-TARE function simply subtracts the added weight of the chair or wheelchair and provides the weight of the patient alone. Thanks to its high load-bearing capacity of 300 kg, the scale is perfectly suited for use by seriously overweight persons. The ramp seca 470 is included in standard delivery.

Versatile wired remote display.


Because the display element of the seca 675 is connected to the scale by means of 2.30 m long cable, it can be positioned where the user can most easily read the results. A wall mounting is also included in the standard delivery.

Optional: non-slip base.

The seca 472 stand is the stable solution for cable remote displays of seca 675. The stable construction with non-slip rubber feet replaces the former model on transport castors and can be positioned everywhere quickly and securely. The display element can easily be attached to the stand by means of a stable bracket. The height of the display is continuously adjustable up to 100 cm.

Easy to transport.

Despite its size, the seca 675 is amazingly easy to handle. When upright, the scale can be moved easily on castors and steered by means of the integrated handle. This feature makes the seca 675 suitable for mobile use.

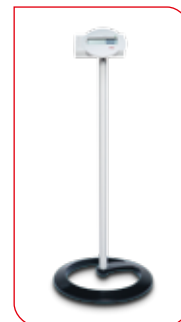

Stand for cable remote displays of seca scales and measuring rods.



A ramp for easy roll-on of wheelchair is included in delivery.

seca 675

Technical Data

- Capacity: 300 kg
- Graduation: 100 g < 200 kg > 200 g
- Dimensions (WxHxD): 904 x 70 x 1,064 mm
- Dimensions/platform (WxHxD): 760 x 55 x 930 mm
- Weight: 27 kg
- Power supply: power adapter
- Functions: TARE, pre-TARE, HOLD, auto-HOLD, BMI, auto-CLEAR, auto weighing range switch-over, damping, automatic switch-off
- Optional: stand for cable remote display seca 472, second ramp seca 470
- Approval class:

Medguard Professional Healthcare

Tel: 01 835 2378

Email: orders@medguard.ie

Fax: 01 969 5009

Web: www.medguard.ie

Unit 6/7, Block 4, Ashbourne Business Park, Ashbourne, Co. Meath, Ireland