



# One-piece Guedel airways

Safer by design



**Airway Management** ▪ Airway Devices

## One-piece Guedel airways

### Safer by design

By using the latest manufacturing technology Intersurgical has produced the first major breakthrough in Guedel airway design for thirty years.


Combining two materials with very different physical characteristics into one moulding, our One-piece design eliminates the separate bite block found in all conventional Guedel airways.

Loose or detached bite blocks constitute a real threat to patient safety and have been the subject of several hazard notices in the past. Our One-piece design overcomes all of the associated risks.

As well as improved safety, the One-piece also incorporates a soft tip reducing the risk of trauma to the patient. We also offer a sterile option, providing a further choice.



▶ Video available at [www.intersurgical.com](http://www.intersurgical.com)

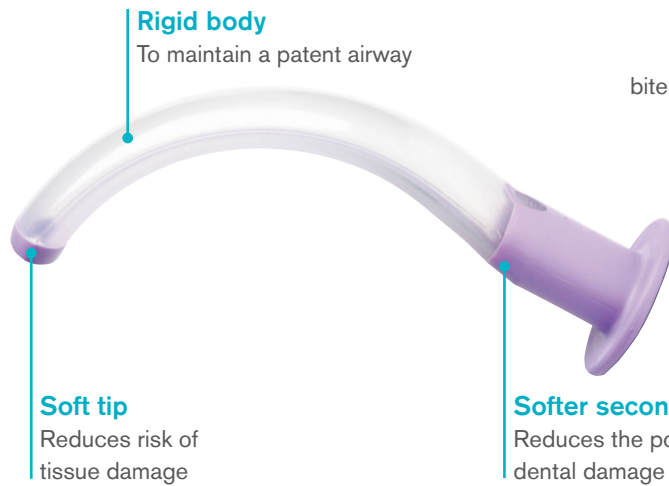










**One-piece design**  
With no separate inner bite block – safer by design

**Nine colour-coded sizes**  
To meet all patient requirements




**Softer second material**  
Reduces the potential for dental damage

**Clear colour-coding**  
For easy identification

Code	Description	Size	ISO	Colour	Box Qty.
1115120 (S*)	One-piece Guedel airway	5	12.0	Purple	70 
1114100 (S*)	One-piece Guedel airway	4	10.0	Red	90 
1113090 (S*)	One-piece Guedel airway	3	9.0	Orange	90 
1112080 (S*)	One-piece Guedel airway	2	8.0	Green	100 
1111570 (S*)	One-piece Guedel airway	1.5	7.0	Yellow	100 
1111065 (S*)	One-piece Guedel airway	1	6.5	White	20 
1110055 (S*)	One-piece Guedel airway	0	5.5	Grey	25 
1100050 (S*)	One-piece Guedel airway	00	5.0	Blue	25 
1000035 (S*)	One-piece Guedel airway	000	3.5	Pink	25 

Make an enquiry

 Lower environmental impact product

 Sterile option available

(S\*) Add an S to the seven digit code number for the sterile version of this product eg. 5000000S



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK

T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 [info@intersurgical.com](mailto:info@intersurgical.com) [www.intersurgical.com](http://www.intersurgical.com)



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

**Please think before you print**

Save energy and paper.  
If you must print this information sheet please print it double sided.

IS6.10 • Issue 3 07.19