



Overview:

Designed for phlebotomy procedures, this fixed height chair features manual synchronised adjustment of the backrest and foot section for patient positioning. This facility enables the clinician to sit the patient comfortably and also provides a flat surface with head-down recovery facility. The chair provides a comfortable sitting position and is supplied with a pair of multi positional phlebotomy arms with height and angle adjustment. The chair is available with fixed non slip feet or braked castors and is ideal for GP Surgeries and health centres etc.




Standard Features

- UK designed & manufactured
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Safe working load 250kg
- Synchronised backrest and foot section
- Fixed height powder coated steel frame
- Adjustable phlebotomy arms
- Hardwearing antimicrobial vinyl
- Choice of fixed feet or twin manual toe-braked castors
- 70mm comfort foam
- Hoist clearance 13cm (5")

Overall Dimensions

Length: 176cm (69") approx
 Width: 63cm (25") approx
 Height: 52cm (20") approx

Upholstery Dimensions

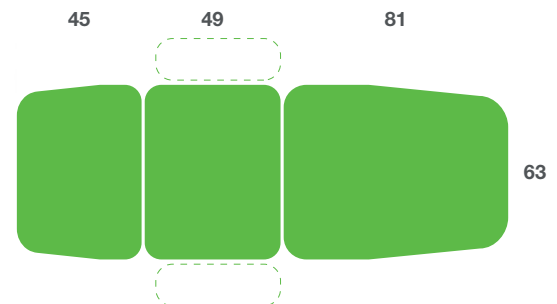

Options & Accessories



Breathing Hole



Paper Roll Holder



D Cushion



IV Pole



Phlebotomy Trolley



Operator's Chair

Operational Parameters



Model 113B - Fixed Feet fitted

Model 113BW - Wheels fitted

Product Benefits

- UK designed & manufactured with long-term support
- Lifetime frame and 5 year parts guarantee providing longevity of product life and value for money
- Safe working load of 250kgs to cater for very heavy patients
- Powder coated fixed height steel frame for long life
- Manual synchronised adjustment of backrest and foot section angles
- The backrest is adjustable from 85° to -15° providing head-down patient recovery
- The foot section is manually adjustable from 10° to -90° to provide patient access from the front or side
- Supplied with a pair of multi-angular adjustable PU arms for blood taking procedures
- The PU arms are easy to clean and feature a sculpted design for patient security and are height adjustable for tall patients
- Deep 70mm comfort foam for increased patient comfort
- Hardwearing, easy to clean antimicrobial upholstery vinyl in a wide choice of colours

Technical Standards

- Safe Working Load 250kg (550lb / 39 Stone)
- IP Rating - IPx4
- EE Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601. P (EN)
- Vinyl flame retardant to CRIB 5 (BS 5852)
- Foam flame retardant to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]
- CE Marked
- Externally audited to ISO9001:2008 by LQRA

Upholstery Vinyl - Phthalate Free

- Flame retardant to CRIB 5
- Independantly tested to verify effectiveness at reducing the spread of MRSA, E.coli and other similar bacteria and fungi
- Rated at 1,000,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high traffic, heavy duty areas



Jasmine



Battleship



Cobalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond